

## CELLBOOSTER® GLOW

### Hyaluronic Acid 18mg Revitalizing Complex:

Non-crosslinked Hyaluronic Acid 18 mg (6 mg/ml) with Amino Acids and Vitamins.



#### DESCRIPTION

CELLBOOSTER® GLOW: Hyaluronic Acid Revitalizing Complex is produced using patented CHAC Technology.

It consists of non-crosslinked Hyaluronic Acid with Amino Acids and Vitamins. The device is designed for nourishment, moisturizing and deep hydration which reduces hyperpigmentation and fights against the oxidative stress.

It restores dermal structural integrity, neutralizes free radicals and reduces the effects of sun damage.

The device is a clear, colorless, sterile, injectable solution for single use presented in a case of 6 × 3 mL vials.



## COMPOSITION

Each vial contains:

- Non-crosslinked Sodium Hyaluronate 18 mg
- Amino acids (Cysteine, Glutathione, Glycine, Lysine, Proline, Valine)
- Vitamins (Sodium Ascorbyl Phosphate, Biotin)
- Phosphate buffer pH 6.8-7.6 q.s.



## INDICATIONS:

- Hyperpigmentation.
- First signs of aging, wrinkles.
- Dry skin.
- Sun-damaged and environmentally compromised skin.
- Uneven skin tone

## FORMAT:

Box with 6 vials x 3 ml.

## PROGRAM:

### Intensive course:



3 treatments with 2-3 weeks interval.

### Support course:



1 treatment per 2 months.

## TREATMENT AREAS:


Face, neck, décolleté area, back of the hands, internal face of the arms or body.

## PROTOCOL A): Full face treatment with CELLBOOSTER® GLOW

### SKIN PREPARATION

- Cleansing | • Antiseptic | • Topical anesthesia (optional)

### 1 PAPULA

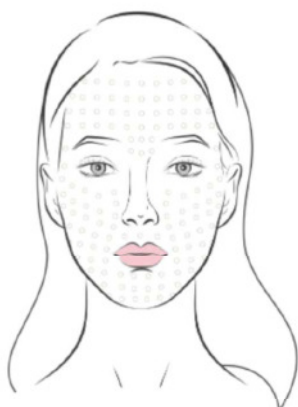

- **Injection angle:** 30°, bevel of the needle up.
- **Injection depth:** Intradermal.
- **Volume injected:** 3ml for full face, 0,025 ml per point every 2/3mm.
- **Needle type:** 30-32G, 4 mm.

### 2 DEEP PAPULA



- **Injection angle:** 45°, bevel of the needle down.
- **Injection depth:** 2mm, on the pigmentation area.
- **Volume injected:** 0,025 ml per injection
- **Needle type:** 30-32G, 4 mm.

### 3 MICROPAPULA



- **Injection angle:** 30°, bevel of the needle up.
- **Injection depth:** Intradermal, Superficial.
- **Volume injected:** 0,025 ml per point under each eye.
- **Needle type:** 32G or smaller, 4 mm.

### 4 NAPPAGE

- **Injection angle:** 30-35°.
- **Injection depth:** epidermal.
- **Volume injected:** minimal.
- **Needle type:** 30G/smaller.



### FINALIZATION


- Antiseptic.
- Post Injection cream.

## PROTOCOL B): Decollete and neck areas treatment with CELLBOOSTER<sup>®</sup> GLOW

### SKIN PREPARATION

•Cleansing | •Antiseptic | •Topical anesthesia (optional)

### 1 PAPULA

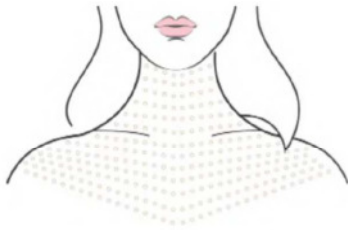

- **Injection angle:** 30°, bevel of the needle up.
- **Injection depth:** intradermal.
- **Volume injected:** 3ml for neck and décolletage, 0,025 ml per point every 2/3 mm.
- **Needle type:** 30-32G, 4 mm.


### 2 DEEP PAPULA



- **Injection angle:** 45°, bevel of the needle down.
- **Injection depth:** 2mm, on pigmentation or redness area.
- **Volume injected:** 0,025 ml per injection
- **Needle type:** 30-32G, 4 mm.

### 3 RETROGRADE LINEAR

- **Injection angle:** 30°, then parallel to the skin.
- **Injection depth:** Intradermal.
- **Volume injected:** 0,02-0,05 ml on each line.
- **Needle type:** 30G or smaller, 13 mm.

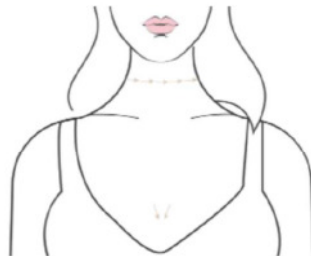

### FINALIZATION

- Antiseptic
- Post injection cream